

DATA VALIDATION REPORT

**Red Hill Bulk Fuel Storage Facility
Joint Base Pearl Harbor-Hickam
CV 23F0104**

**SDG: FC11160
SGS North America Inc - Orlando**

Prepared by
ENVIRONMENTAL DATA SERVICES, LTD.


Prepared for
AECOM Environmental

Released: 12/01/23


Data Validators and Peer Reviewers:




Diane Waldschmidt

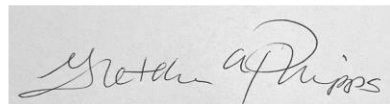

Gretchen Phipps

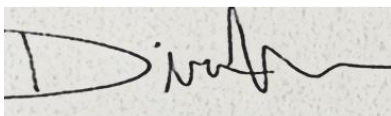

Dina Manov

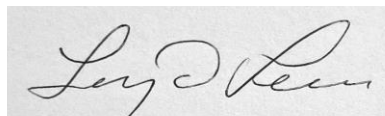

Larry Lewis



Paloma Hoelzle

EXECUTIVE NARRATIVE

Sample Delivery Group: FC11160

Laboratory: SGS North America Inc - Orlando

Site: Red Hill Bulk Storage Facility, CV 23F0104

Sampling dates: 11/09/2023

Number of Samples: 2

Test Method: USEPA Method 1633

Analysis: per- and polyfluoroalkyl substances (PFAS)

Quality Assurance Project Plan: Sampling and Analysis Plan, Investigation and Remediation of Releases and Groundwater Protection and Evaluation, Red Hill Bulk Fuel Storage Facility, Joint Base Pearl Harbor-Hickam, O'ahu, Hawai'i (Revision 01, April 2017); PFAS-Specific Sampling and Analysis plan, Red Hill Bulk Fuel Storage Facility, Adit 6, Joint Base Pearl Harbor-Hickam, O'Ahū, Hawai'i (November 30, 2022) (SAP).

Validation Guidelines: United States Department of Defense Data Validation Guidelines Module 6: Data Validation Procedure for Per- and Polyfluoroalkyl Substances analysis by QSM Table B-24, Environmental Data Quality Workgroup, October 18, 2022; United States Department of Defense (DOD) Environmental Data Quality Workgroup (EDQW), General Validation Guidelines, November 2019.

| Client Sample Identification | Laboratory Sample Identification | Matrix | Validation Stage |
|------------------------------|----------------------------------|-------------|------------------|
| AF-RHMW12A-WGN01LF-2311 | FC11160-1 | groundwater | S2BVEM |
| AF-RHMW12A-WGFD01LF-2311 | FC11160-2 | groundwater | S2BVEM |

Table 1 provides a summary of the major and minor data quality issues identified in this data set. All data are acceptable except those results which have been qualified with "X", rejected. Data validation qualifiers along with associated descriptions are provided in Table 2. All data qualification related to this group of samples is detailed on the attached sheets.

All data users should note two facts. First, an "X" flag means that the associated value is unusable due to significant quality control (QC) problems, the data is invalid and provides no information as to whether the compound is present or not. "X" values should not appear on any data tables even as a last resort. Second, no analyte concentration, even if it passed all QC tests, is guaranteed to be accurate. Strict QC serves to increase confidence in data, but any value potentially contains error.

DATA ASSESSMENT

1. NARRATIVE AND COMPLETENESS REVIEW

The case narrative was reviewed, and the data package was checked for completeness. No discrepancies were noted.

2. SAMPLE DELIVERY AND CONDITION

The samples arrived at the laboratory in acceptable condition. Proper custody was documented.

3. HOLDING TIME

The amount of an analyte in a sample can change with time due to chemical instability, degradation, volatilization, etc. If the specified holding time is exceeded, the data may not be valid. Proper sample handling and preservation also play a role in the chemical stability of analytes in the sample matrix. If samples are not collected and stored using proper containers and/or preservatives, data may not be valid.

No problems were found for this criterion.

4. CALIBRATION

Satisfactory instrument calibration is established to ensure that the instrument can produce acceptable quantitative data. An initial calibration demonstrates that the instrument can give acceptable performance at the beginning of an experimental sequence. The continuing calibration checks document that the instrument is giving satisfactory daily performance. Additionally, a continuing calibration is analyzed at the end of each 12-hour analytical sequence, denoted as a "closing" calibration verification and ascertains acceptable performance at the conclusion of the analytical sequence.

A) Initial Calibration

Percent relative standard deviation (%RSD) is calculated from the initial calibration and is used to indicate stability of a specific compound over the calibration range.

An RSD value outside the initial calibration limit indicates the potential for quantitation errors. For this reason, all positive and non-detected results are qualified as estimated. Severe performance failures (RSD >30%) requires rejection of all results. The following QC criteria have been applied for this project: The %RSD of initial calibration must be <20%.

No problems were found for this criterion.

B) Continuing Calibration

The Percent Recovery (%R) for all target analytes in the continuing calibration must be within 70-130%. All initial calibration verification (ICV) and continuing calibration verification (CCV) %Rs were with acceptance limits with the following exceptions.

No problems were found for this criterion.

C) Instrument Sensitivity Check

Prior to analysis an instrument sensitivity check (ISC) must be performed. The ISC must be at the limit of quantitation (LOQ). All analyte concentrations must be within $\pm 30\%$. Note: the laboratory reports refer to the ISC as Low-Concentration Calibration Verification (LCCV). The validator has determined that the low level CCV in the laboratory's report is equivalent to the method required ISC.

No problems were found for this criterion.

5. BLANK CONTAMINATION

Quality assurance (QA) blanks, i.e., method, field, or rinse blanks are prepared to identify any contamination which may have been introduced into the samples during sample preparation or field activity. Method blanks measure laboratory contamination. Field and rinse blanks measure cross-contamination of samples during field operations. When an equipment blank, or lab blank has an analyte detection, then all associated field samples are qualified per validation guidance as appropriate.

A) Method blank contamination:

No problems were found for this criterion.

B) Instrument blank contamination:

No problems requiring result qualification were found for this criterion.

C) Field/Equipment blank contamination:

No samples were submitted as a field/equipment blank in association with samples in this sample delivery group (SDG).

6. EXTRACTED INTERNAL STANDARDS

All samples are spiked with labeled standard compounds prior to sample preparation and analyses to evaluate overall laboratory performance and efficiency of the analytical technique. The reported project samples observed surrogate recoveries within the established limits in all cases with the following exceptions.

No problems were found for this criterion.

7. NON-EXTRACTED INTERNAL STANDARDS

Non-extracted internal standard peak areas are used to quantify extracted internal standard recoveries. The reported project samples had non-extracted internal standard area counts within the established limits in all cases with the following exceptions.

No problems were found for this criterion.

8. COMPOUND IDENTIFICATION

The project target analyte compounds are identified on the LC/MS/MS by using the analytes retention time (RT). The retention time of each target analyte should be within ± 0.4 minutes of the predicted retention. Target analyte detections should display a signal-to-noise of $\geq 3:1$, have proper peak integration, and display all ions at the correct retention times.

Target analyte detections should have passing ion ratios (50 - 150% of theoretical). Ion ratio failures could be caused by matrix interference and/or be the result of the presence of isomers in the sample at different ratios than the ratio of isomers present in the calibration standards.

Target compound identification was verified. No anomalies were identified.

9. COMPOUND QUANTIFICATION

Target compound quantitation was not verified as part of the Level 2B data validation.

Manual integrations were not reviewed at the Stage 2B level.

10. MATRIX SPIKE/MATRIX SPIKE DUPLICATE RECOVERY / MATRIX DUPLICATE

Matrix spike/matrix spike duplicate (MS/MSD) data are generated to determine the long-term precision and accuracy of the analytical method in various matrices. The MS/MSD data may be used in conjunction with other quality control criteria for additional qualification of data.

No sample was submitted for MS/MSD and/or matrix duplicate evaluation in association with this SDG.

11. FIELD DUPLICATES/TRIPPLICATES

Field duplicates may be taken and analyzed as an indication of overall precision. These analyses measure both field and laboratory precision. A control limit of $\leq 30\%$ for the Relative Percent Difference (RPD) for water samples and $\leq 50\%$ RPD for solid samples, shall be used when original and duplicate sample values are greater than or equal to the sample specific LOQ. For field duplicate analyses that do not meet the technical criteria, the action was applied to only the parent sample and its duplicate. A control limit of $\leq 35\%$ RSD was applied for field triplicate samples when original and triplicate sample values are greater than the sample specific LOQ. For field triplicate analyses that do not meet the technical criteria, the action was applied to only the parent sample, duplicate and triplicate.

Samples AF-RHMW12A-WGN01LF-2311 and AF-RHMW12A-WGFD01LF-2311 were submitted as a field duplicate pair in association with this SDG. Upon evaluation adequate precision was demonstrated.

12. LABORATORY CONTROL SAMPLES

The Laboratory Control Sample (LCS) serves as a monitor of the overall performance of each step during the analysis, including the sample preparation. The LCS results are used to verify that the laboratory can perform the analysis in a clean matrix. Note: in addition to the standard LCS the laboratory has also provided a second LCS referred to as the MRL check in the laboratory report.

No problems requiring result qualification were found for this criterion.

13. DILUTIONS, RE-EXTRACTIONS & REANALYSIS

Samples may be re-analyzed for dilution, re-extraction and for other QC reasons. In such cases, the best result values are used.

No sample re-extracted and/or reanalyzed were provided for data review.

14. SYSTEM PERFORMANCE AND OVERALL ASSESSMENT

Overall, the laboratory data generated met the project goals and quality control criteria, with the exceptions identified in this report and as summarized in Table 1.

**Table 1
Review Elements Summary**

| | Were acceptance criteria met? | | |
|--|-------------------------------|-------|-------|
| | Yes | No | |
| | | Major | Minor |
| Per-fluorinated Compounds | | | |
| Holding Time/Sample Handling | x | | |
| Method Blanks | x | | |
| Instrument Blanks | x | | |
| Field Blanks | NA | | |
| Calibration Percent Relative Standard Deviation and Percent Difference | x | | |
| Instrument Sensitivity Check | x | | |
| Extracted Internal Standards | x | | |
| Non-Extracted Internal Standards | x | | |
| Compound Identification | x | | |
| Matrix Spike/Matrix Spike Duplicate/Matrix Duplicate | NA | | |
| Laboratory Control Sample | x | | |
| Other Quality Control Data out of Specification | x | | |
| Field Duplicate / Triplicate | x | | |

Major= Major data quality issue identified resulting in rejection of data.

Minor= Minor data quality issue identified resulting in the qualification of data. Data qualification should be used to inform the data users of data limitations.

NA = Not applicable

**Table 2
Data Validation Qualifiers**

| Data Qualifier | Definition |
|-----------------------|---|
| U | The analyte was analyzed for but was not detected above the level of the reported sample quantitation limit. |
| J | The result is an estimated quantity. The associated numerical value is the approximate concentration of the analyte in the sample. |
| J+ | The result is an estimated quantity, but the result may be biased high. |
| J- | The result is an estimated quantity, but the result may be biased low. |
| UJ | The analyte was analyzed for but was not detected. The reported quantitation limit is approximate and may be inaccurate or imprecise. |
| X | The sample results (including non-detects) were affected by serious deficiencies in the ability to analyze the sample and to meet published method and project quality control criteria. The presence or absence of the analyte cannot be substantiated by the data provided. |
| R | The data are unusable. The sample results are rejected due to serious deficiencies in meeting QC criteria. The analyte may or may not be present in the sample. |

**Table 3
PFAS Definitions Table**

| NO | CAS # | Target Name | Target Abbreviation |
|-----------|--------------|---|----------------------------|
| 1 | 763051-92-9 | 11-Chloroeicosafluoro-3-oxaundecane-1-sulfonic acid | 11Cl-PF3OUdS |
| 2 | 914637-49-3 | 2H,2H,3H,3H-Perfluorooctanoic acid | 5:3FTCA |
| 3 | 812-70-4 | 3-Perfluoroheptyl propanoic acid | 7:3FTCA |
| 4 | 356-02-5 | 3-Perfluoropropyl propanoic acid | 3:3FTCA |
| 5 | 919005-14-4 | 4,8-Dioxa-3H-perfluorononanoic acid | ADONA |
| 6 | 757124-72-4 | 4:2 Fluorotelomer sulfonic acid | 4:2 FTS |
| 7 | 27619-97-2 | 6:2 Fluorotelomer sulfonic acid | 6:2 FTS |
| 8 | 39108-34-4 | 8:2 Fluorotelomer sulfonic acid | 8:2 FTS |
| 9 | 756426-58-1 | 9-Chlorohexadecafluoro-3-oxanone-1-sulfonic acid | 9Cl-PF3ONS |
| 10 | 13252-13-6 | Hexafluoropropylene oxide dimer acid | HFPO-DA |
| 11 | 4151-50-2 | N-Ethyl perfluorooctanesulfonamide | NEtFOSA |
| 12 | 2991-50-6 | N-Ethyl perfluorooctanesulfonamidoacetic acid | NEtFOSAA |
| 13 | 1691-99-2 | N-Ethyl perfluorooctanesulfonamidoethanol | NEtFOSE |
| 14 | 31506-32-8 | N-Methyl heptadecafluorooctanesulfonamide | NMeFOSA |
| 15 | 2355-31-9 | N-Methyl perfluorooctanesulfonamidoacetic acid | NMeFOSAA |
| 16 | 24448-09-7 | N-Methyl perfluorooctanesulfonamidoethanol | NMeFOSE |
| 17 | 151772-58-6 | Nonafluoro-3,6-dioxaheptanoic acid | NFDHA |
| 18 | 113507-82-7 | Perfluoro(2-ethoxyethane)sulfonic acid | PFEESA |
| 19 | 377-73-1 | Perfluoro-3-methoxypropanoic acid | PFMPA |
| 20 | 863090-89-5 | Perfluoro-4-methoxybutanoic acid | PFMBA |
| 21 | 375-73-5 | Perfluorobutanesulfonic acid | PFBASA |
| 22 | 375-22-4 | Perfluorobutanoic acid | PFBA |
| 23 | 335-77-3 | Perfluorodecanesulfonic acid | PFDS |
| 24 | 335-76-2 | Perfluorodecanoic acid | PFDA |
| 25 | 79780-39-5 | Perfluorododecanesulfonic acid | PFDoS |
| 26 | 307-55-1 | Perfluorododecanoic acid | PFDoA |
| 27 | 375-92-8 | Perfluoroheptanesulfonic acid | PFHpS |
| 28 | 375-85-9 | Perfluoroheptanoic acid | PFHpA |
| 29 | 355-46-4 | Perfluorohexanesulfonic acid | PFHXSA |
| 30 | 307-24-4 | Perfluorohexanoic acid | PFHxA |
| 31 | 68259-12-1 | Perfluorononanesulfonic acid | PFNS |
| 32 | 375-95-1 | Perfluorononanoic acid | PFNA |
| 33 | 754-91-6 | Perfluorooctanesulfonamide | PFOSA |
| 34 | 1763-23-1 | Perfluorooctanesulfonic acid | PFOS |
| 35 | 335-67-1 | Perfluorooctanoic acid | PFOA |
| 36 | 2706-91-4 | Perfluoropentanesulfonic acid | PFPeS |
| 37 | 2706-90-3 | Perfluoropentanoic acid | PFPeA |
| 38 | 376-06-7 | Perfluorotetradecanoic acid | PFTeDA |
| 39 | 72629-94-8 | Perfluorotridecanoic acid | PFTTrDA |
| 40 | 2058-94-8 | Perfluoroundecanoic acid | PFUnA |

| Data Qualification Reason Codes | |
|--|--|
| Reason Code | Reason Code Description |
| A | Serial dilution |
| A1 | Ambient Blank |
| B | The analyte was found in an associated blank as well as in the sample. |
| B2 | CCB |
| B3 | CCB - Neg |
| B4 | Grinding Blank |
| C | LCS Recovery |
| C1 | Reference Recovery |
| C2 | Reference Recovery RPD |
| D | MS RPD |
| D1 | Lab Replicate RPD |
| D2 | No precision available |
| D3 | Field Duplicate RPD |
| D4 | Field Triplicate RSD |
| D5 | Laboratory Triplicate RSD |
| F | Field Blank |
| F1 | Hydrocarbon pattern does not match standard |
| G1 | Initial Calibration RRF |
| G2 | Initial Calibration RSD/r²/r |
| G3 | ICV RRF |
| H1 | Test Hold Time |
| H2 | Prep Hold Time |
| I | Surrogate recovery outside project limits. |
| J | CRA/CRI Recovery |
| K | An analyte (non-common laboratory artifact) was detected in the sample at a concentration less than 5X the concentration detected in the associated method blank. |
| L | Lab Blank |
| L1 | Lab Blank - Neg |
| M | MS Recovery |
| M2 | Post Spike |
| N | Blank - No Action |
| O | ICS |
| P | Sample preservation/collection requirement not met. |
| P1 | Column RPD |
| P2 | Improper preparation/extraction |
| Q | Encore sample holding time exceeded by more than 2X. |
| Q1 | Material Blank |

| | |
|----|--|
| Q2 | Encore sample holding time exceeded by less than 2X. |
| R | Exceeds Linear Calibration Range |
| S | Internal standard |
| T | Trip Blank |
| TI | Tentatively Identified Compound |
| TR | Trace Level Detect |
| U | Receipt Temperature |
| V | Equipment Blank |
| V1 | ICV |
| V2 | CCV |
| V3 | CCV RRF |
| V4 | Sample Receipt Condition |
| V5 | Ending Continuing Calibration Verification |
| V6 | Low Level Calibration Verification |
| V7 | Interference Check Sample A |
| V8 | Interference Check Sample AB |
| V9 | Interference Check Sample A - Negative |
| W | Column breakdown (pesticides/8270) |
| X | Raised reporting limit |
| Y | Cooler temperature greater than 10 degreeec C. |
| Y1 | False Positive |
| Y2 | Data rejected due to radiological anomolies |
| Y3 | Non-accredited analyte/compound. Accreditation not offered at time of analyses for the analyte/compound by the stated method and matrix. |
| Y4 | Performance Check - Degradation of DDT |
| Y5 | Extracted Internal Standard |
| Y6 | Analyte not confirmed on second column. |
| Y7 | Signal to Noise Ratio not met |
| Z | LCS RPD |
| Z1 | Non-accredited analyte/compound |
| Z1 | Data rejected, more valid data available. |
| Z2 | Detection Level not met uncertainty greater than DL |
| Z4 | MDA Greater than RDL. |
| Z5 | Ion Ratio |
| Z6 | Samples were analyzed past the 12 hour time period from the Tune or opening CCV. |

Calculation Documentation

Field Duplicate Calculations
 RPD Water <30%
 RPD Solid <50%

| Parameters | Original Sample | | Duplicate Sample | | RPD | LOQ |
|--|-------------------------|--|--------------------------|--|---------|-----|
| | AF-RHMW12A-WGN01LF-2311 | | AF-RHMW12A-WGFD01LF-2311 | | | |
| 11-Chloroicosafuoro-3-oxaundecane-1-sulfonic acid (11Cl-PF3OUdS) | 0 | | 0 | | #DIV/0! | |
| 2H,2H,3H,3H-Perfluorooctanoic acid (5:3FTCA) | 0 | | 0 | | #DIV/0! | |
| 3-Perfluorohexyl propanoic acid (7:3FTCA) | 0 | | 0 | | #DIV/0! | |
| 3-Perfluoropropyl propanoic acid (3:3FTCA) | 0 | | 0 | | #DIV/0! | |
| 4,8-Dioxa-3H-perfluorononanoic acid (ADONA) | 0 | | 0 | | #DIV/0! | |
| 4:2 Fluorotelomer sulfonic acid (4:2 FTS) | 0 | | 0 | | #DIV/0! | |
| 6:2 Fluorotelomer sulfonic acid (6:2 FTS) | 0 | | 0 | | #DIV/0! | |
| 8:2 Fluorotelomer sulfonic acid (8:2 FTS) | 0 | | 0 | | #DIV/0! | |
| 9-Chlorohexadecafluoro-3-oxanone-1-sulfonic acid (9Cl-PF3ONS) | 0 | | 0 | | #DIV/0! | |
| Hexafluoropropylene oxide dimer acid (HFPO-DA) | 0 | | 0 | | #DIV/0! | |
| N-Ethyl perfluorooctanesulfonamide (NEtFOSA) | 0 | | 0 | | #DIV/0! | |
| N-Ethyl perfluorooctanesulfonamidoacetic acid (NEtFOSAA) | 0 | | 0 | | #DIV/0! | |
| N-Ethyl perfluorooctanesulfonamidoethanol (NEtFOSE) | 0 | | 0 | | #DIV/0! | |
| N-Methyl heptadecafluorooctanesulfonamide (NMeFOSA) | 0 | | 0 | | #DIV/0! | |
| N-Methyl perfluorooctanesulfonamidoacetic acid (NMeFOSAA) | 0 | | 0 | | #DIV/0! | |
| N-Methyl perfluorooctanesulfonamidoethanol (NMeFOSE) | 0 | | 0 | | #DIV/0! | |
| Nonafluoro-3,6-dioxaheptanoic acid (NFDHA) | 0 | | 0 | | #DIV/0! | |
| Perfluoro(2-ethoxyethane)sulfonic acid (PFEEESA) | 0 | | 0 | | #DIV/0! | |
| Perfluoro-3-methoxypropanoic acid (PFMPA) | 0 | | 0 | | #DIV/0! | |
| Perfluoro-4-methoxybutanoic acid (PFMBA) | 0 | | 0 | | #DIV/0! | |
| Perfluorobutanesulfonic acid (PFBS) | 0 | | 0 | | #DIV/0! | |
| Perfluorobutanoic acid (PFBA) | 0 | | 0 | | #DIV/0! | |
| Perfluorodecanesulfonic acid (PFDS) | 0 | | 0 | | #DIV/0! | |
| Perfluorodecanoic acid (PFDA) | 0 | | 0 | | #DIV/0! | |
| Perfluorododecanesulfonic acid (PFDoS) | 0 | | 0 | | #DIV/0! | |
| Perfluorododecanoic acid (PFDoA) | 0 | | 0 | | #DIV/0! | |
| Perfluoroheptanesulfonic acid (PFHpS) | 0 | | 0 | | #DIV/0! | |
| Perfluoroheptanoic acid (PFHpA) | 0 | | 0 | | #DIV/0! | |
| Perfluorohexanesulfonic acid (PFHxS) | 0 | | 0 | | #DIV/0! | |
| Perfluorohexanoic acid (PFHxA) | 0 | | 0 | | #DIV/0! | |
| Perfluoronanesulfonic acid (PFNS) | 0 | | 0 | | #DIV/0! | |
| Perfluorononanoic acid (PFNA) | 0 | | 0 | | #DIV/0! | |
| Perfluorooctanesulfonamide (PFOSA) | 0 | | 0 | | #DIV/0! | |
| Perfluorooctanesulfonic acid (PFOS) | 0 | | 0 | | #DIV/0! | |
| Perfluorooctanoic acid (PFOA) | 0 | | 0 | | #DIV/0! | |
| Perfluoropentanesulfonic acid (PFPeS) | 0 | | 0 | | #DIV/0! | |
| Perfluoropentanoic acid (PFPeA) | 0 | | 0 | | #DIV/0! | |
| Perfluorotetradecanoic acid (PFTeDA) | 0 | | 0 | | #DIV/0! | |
| Perfluorotridecanoic acid (PFTriDA) | 0 | | 0 | | #DIV/0! | |
| Perfluoroundecanoic acid (PFUnA) | 0 | | 0 | | #DIV/0! | |

Flag

Data Validation Worksheet

DATA VALIDATION PFAS

Module 6; PFAS by QSM Table B-24; October 18, 2022

Validator: [GAP](#)

Reviewer: [DLW](#)

Date Validated: [12/1/23](#)

Reviewed: [12/1/23](#)

Project: Red Hill

SDG: FC11160

LAB: SGS North America Inc. - Orlando

Samples Collected: 11/9/2023

[2 GW](#)

SAMPLE RECEIPT AND CASE NARRATIVE REVIEW

- ✓ Traffic reports, chain-of-custody forms or SDG narrative do not indicate any problems with sample receipt, condition of the samples, analytical problems or special circumstances affecting the quality of the data.
- ✓ AFFF samples are to be shipped in HDPE containers with an unlined cap
- ✓ Shipment temp 0-6°C: recommended to freeze tissue samples upon receipt
- ✓ If temp upon receipt is greater than 6°C J/UJ all

[Received on 11/10 at 5 C](#)

HOLDING TIMES

- ✓ Recommended storage temp is $\leq -20^{\circ}\text{C}$
- ✓ Per method 1633: aqueous samples may be held in the lab for up to 90 days when stored at recommended temp and protected from light; when stored at 0-6 °C and protected from light samples can be held for up to 28 days (see method for additional details)
- ✓ Per method 1633: solid samples may be held in the lab for up to 90 days when stored at recommended temp or 0-6 °C (see method for additional details)
- ✓ Per method 1633: biosolid samples may be held in the lab for up to 90 days when stored at recommended temp or 0-6 °C; however, freezing is recommended (see method for additional details)
- ✓ Samples extracts should be stored at 0-4°C protected from light and analyzed within 90 days
- ✓ If hold time is exceeded qualify J/UJ

- ✓ If hold time is grossly exceeded (2X hold time) J/X

244 **Table II. Sample Storage and Holding Time Requirements**

| Matrix Type | Stored at 0 - 6°C, protected from light | | Stored at ≤ -20°C, protected from light | |
|------------------|---|--|---|---|
| | Holding Time | Caveat | Holding Time | Caveat |
| Aqueous | 28 days | Precursor degradation occurs after 7 days | 90 days | None |
| Solid and Tissue | 90 days | Should be prepared as soon as possible if NFDHA is a target analyte | 90 days | Should be prepared as soon as possible if NFDHA is a target analyte |
| Biosolid | 90 days | Not recommended due to the production of gases due to microbiological activity | 90 days | None |

Samples collected 11/09/2023

Extracted 11/13

Analyzed 11/15

All ok

Extracted Internal STANDARDS

- ✓ Added to all QC and field samples
- ✓ Recoveries are within the limits as defined in QAPP; otherwise QSM criteria (20-150%) should be used
- ✓ Detected for analytes qualified using an EIS percent recovery >200% should be qualified J-. Nod-detects should not be qualified.
- ✓ If EIS recovery is <10%; associated detected and non-detects should be qualified X
- ✓ EIS retention times should be within 0.4 minutes of standard; use professional judgment to qualify

For Red Hill project(see Kristin's email on file in project folder 12/14/22 at 3:25pm)

For EIS %Rs >150% J- positive results, no action on non-detects

For EIS %Rs between lab limit of 20-150%; no action

For EIS %Rs <20% but >10%; J+ positive results, UJ non-detects

For EIS %Rs <10% X positive and non-detected (and recommend R of non-detected, J+ of positive results)

ALL OK

Non-Extracted Internal STANDARDS

- ✓ Used to quantify EIS
- ✓ If low are counts are reported (<30%) detected and non-detected should be qualified X

All ok

Laboratory Control Sample (LCS) and Low-Level Laboratory Control Sample (LLCS)

- ✓ LCMS Lab Control Recovery (Form III), Form I, prep log, run log
- ✓ LCS prepared, extracted, analyzed, and reported once for every 20 field samples of a similar matrix, per SDG.
- ✓ Laboratory Control Samples were analyzed for all the target analytes that the samples are analyzed for.
- ✓ Use limits as defined in QAPP; otherwise lab limits or QSM criteria of 40-150%.
- ✓ If LCS or LLCS %R is > upper limit; qualify detects J+; no action on non-detected
- ✓ If LCS or LLCS %R is < lower limit; qualify detected J- and non-detected X

Use lab limits (40-150) to evaluate
All 40 compounds included.

OP58-LLBS all ok

OP58-BS associated with all samples
3:3 Fluorotelomer carboxylate ↑ all samples ND; no action

MS/MSD and Matrix Duplicate

- ✓ LCMS Matrix Spike Recovery (Form III)
- ✓ The Matrix Spike Samples were spiked and analyzed for all the target analytes that the samples are analyzed for (Same analytes as LCS).

- ✓ Per module 6: MS and MSD are applicable where the spike concentration is a least 3 times greater than the native analyte concentration (**3X rule**)
- ✓ Use limits as defined in QAPP; otherwise lab limits or QSM criteria of 40-150%.
- ✓ If MS or MSD %R is > upper limit; qualify detects J+; no action on non-detected
- ✓ If MS or MSD %R is < lower limit but >10%; qualify detected J- and non-detected UJ
- ✓ If MS or MSD %R is < 10%; qualify detected J- and non-detected X
- ✓ If MS/MSD RPD is out; qualify detected J and non-detected UJ
- ✓ For matrix duplicate; for concentrations of analytes that are equal to or greater than the LOQ, the RPD must be ≤30%; if out qualified detected J; no action on non-detects

Use lab limits to evaluate

Sample:

none

Matrix duplicate:

none

BLANKS

- ✓ LCMS Method Blank Summary (Form IV), method blank Form I, prep log, run log
- ✓ Frequency of Analysis: method blank has been analyzed for every 20 (or less) samples of similar matrix or concentration or each extraction batch.
- ✓ Continuing Calibration Blanks (Form I) and run log
- ✓ Frequency of Analysis: immediately following the highest standard analyzed and daily prior to sample analysis.
- ✓ Field/rinse blanks are non-detected for all analytes

312 **Table III: Sample Qualification in the Presence of Blank Contamination**

| | Sample | | |
|------------|----------------------------|-------------------------|----------------------|
| Row Number | Result | Validated Result | Validation Qualifier |
| 1 | Non-detect or detect ≤ LOD | Report at LOD | U |
| 2 | > LOQ but ≤ 5x blank | Report at Sample Result | J+ |
| 3 | > LOQ and > 5x blank | Report at Sample Result | None |

313 LOD = Limit of Detection

All instrument blanks: see calibration section

EB/FBs

none

MASS CALIBRATION

- ✓ Verified to be ± 0.2 amu of true value

Bile Salt Interference Check and Qualitative Identification Standard

- ✓ Provided and requirements met
- ✓ See Module 6

All acceptable

ICAL

- ✓ Initial Calibration Data Curve Evaluation (Form VI) and run log
- ✓ Lowest standard should be at or below LOQ
- ✓ %RSD <20% or relative standard error (RSE) <20%
- ✓ If %RSD > 20% but <30% J/UJ
- ✓ If %RSD >30% J/R

See below

INSTRUMENT PERFORMANCE CHECK PER DRAFT METHOD 1633

- ✓ Concentration equal to LOQ
- ✓ Analyzed after ICAL and daily before samples
- ✓ If not analyzed all associated data should be qualified X
- ✓ The %R for ICV and CCV 30%; if out >130% qualify positive J+ and nondetected UJ; if out <70% qualify positives J- and nondetects UJ
- ✓ Per module if gross exceedances of recoveries <50% or >150%; qualify all associate data X

CCAL

- ✓ Continuing Calibration Data (Form VII) and run log
- ✓ Continuing calibration standard analyzed on each working day, prior to sample analyses.
- ✓ Calibration verification/continuing calibration standard been analyzed after every 10 samples and at the end of each analytical sequence

- ✓ If not analyzed all associated data should be qualified X
- ✓ The %R for ICV and CCV 30%; if out >130% qualify positive J+ and nondetected UJ; if out <70% qualify positives J- and nondetects UJ
- ✓ Per module if gross exceedances of recoveries <50% or >150%; qualify all associate data X

1.0LL CCV is the method required ISC

Instrument GCMS4Q

11/13/2023

all %RSD <20%

S4Q785-ICV785 4Q53740.D 11/13/23 18:12 n/a Initial cal verification 4
 S4Q785-ICV785 4Q53741.D 11/13/23 18:27 n/a Initial cal verification 20

S4Q786-IBLK 4Q53867.D 11/15/23 10:45 n/a Instrument Blank
 4151-50-2 EtFOSA 0.0018 J ug/L ND in all samples; no action
 24448-09-7 MeFOSE 0.0049 J ND in all samples; no action
 S4Q786-CC785 4Q53868.D 11/15/23 11:00 n/a Continuing cal 4
 S4Q786-CC785 4Q53869.D 11/15/23 11:14 n/a Continuing cal 1.0LL

S4Q786-CC785 4Q53879.D 11/15/23 13:53 n/a Continuing cal 4
 S4Q786-ICCB 4Q53880.D 11/15/23 14:08 n/a Continuing Calibration Blank

All samples

S4Q786-CC785 4Q53890.D 11/15/23 16:36 n/a Continuing cal 4

All ok except if noted

COMPOUND IDENTIFICATION

- ✓ RT within ± 0.4 RRT units (review for Level 4)
- ✓ S/N ration 3:1 (review for Level 4)
- ✓ Ion response ratio with $\pm 50\%$ (review for Level 2B)
- ✓ If ion ratio is outside limit; qualify J

Use J flag for module 6

Reason Code: Z5

All ok

FIELD DUPLICATES

- ✓ Use QAPP defined criteria
- ✓ If outside acceptance criteria qualify J/UJ (MODULE FLAGS NONDETECTS TOO)

For field triplicates use 35% RSD per Kristin's email on file from 12/14/22

AF-RHMW12A-WGN01LF-2311 and AF-RHMW12A-WGFD01LF-2311 ok

FIELD DUPLICATE WORKSHEET

Data Validation Report for FC11160

Facility: RH PFAS Sampling
Event: AFFF Assessment Sampling GW 2023 November
SDG: FC11160
Guidance Document: RHS PFAS UFP-QAPP
Prime Contractor: AECOM, Honolulu, HI
Project Manager:
Contract Laboratory(ies): SGS North America, Inc., Orlando, FL
Data Review Contractor:
Data Review Level:
Primary Data Reviewer: ,
Date Submitted:

| Field Sample ID | Lab Sample ID | Matrix | Type/Type Code | E1633DR |
|--------------------------|---------------|--------|--------------------|---------|
| AF-RHMW12A-WGFD01LF-2311 | FC11160-2 | Water | Field Duplicate/FD | X |
| AF-RHMW12A-WGN01LF-2311 | FC11160-1 | Water | Field Sample/N | X |

Data Validation Report for FC11160

This report assesses the analytical data quality associated with the analyses listed on the preceding cover page at data validation level. This assessment has been made through a combination of automated data review (ADR) and supplemental manual review, the details of which are described below. The approach taken in the review of this data set is consistent with the requirements contained in the RHS PFAS UFP-QAPP and the additional guidance documents incorporated by reference to the extent possible. Where definitive guidance is not provided, results have been evaluated in a conservative manner using professional judgment.

Sample collection was managed and directed by AECOM, Honolulu, HI; analyses were performed by SGS North America, Inc., Orlando, FL and were reported under sample delivery group (SDG) FC11160. Data have been evaluated electronically based on electronic data deliverables (EDDs) provided by the laboratory, and hard copy data summary forms have also been reviewed during this effort and compared to the automated review output by the reviewers whose signatures appear on the following page. Findings based on the automated data submission and manual data verification processes are detailed in the ADR narrative and throughout this report.

All quality control (QC) elements associated with this SDG have been reviewed by a project chemist in accordance with the requirements defined for the project. This review is documented in the attached Data Review Checklists. The QC elements listed below were supported by the electronic deliverable and were evaluated using ADR processes.

- Extracted Internal Standard
- Field Duplicate RPD
- Lab Blank
- LCS Recovery
- Prep Hold Time
- Test Hold Time

Results of the ADR process were subsequently reviewed and updated as applicable by the data review chemists identified on the signature page. Quality control elements that were not included in the electronic deliverable were reviewed manually and findings are documented within this report. Summaries of findings and associated qualified results are documented throughout this report.

A total of 0 results (0.00%) out of the 80 results (sample and field QC samples) reported are qualified based on review and 0 results (0.00%) have been rejected or deemed a serious deficiency (X qualifier). Trace values, defined as results that are qualified as estimated because they fall between the detection limit and the reporting limit/limit of quantitation, are not counted as qualified results in the above count. The qualified results are detailed throughout this report and discussed in the narrative below, where appropriate.

Data Validation Report for FC11160

Narrative Comments

| Analytical Method | Data Reviewer Comment |
|-------------------|-----------------------|
|-------------------|-----------------------|

Reviewed by , ,

As the Reviewer, I certify that I have performed a data review process in accordance with the requirements of the project guidance document, and have compared the electronic data to the laboratory's hard copy report and have verified the consistency of the reported sample results and method quality control data between the two deliverables.

Data Validation Report for FC11160

Quality Control Outliers for test method E1633DR, Extracted Internal Standard

Method performance for individual samples is demonstrated through spiking activities. All samples are spiked with internal standards compounds prior to sample preparation (EIS). The sample itself may produce effects due to such factors as interferences and high concentrations of analytes. Summary forms were evaluated and compared to electronic data deliverables. EIS results that were outside of the acceptance criteria are listed below.

| Sample ID/ Lab Sample ID | Analyte | Result | Warning Limits | Control Limits | Units | Qualifier | Reason Code | Comment |
|-----------------------------|---|--------|-------------------|-------------------|---------|-----------|----------------|-----------------------|
| LABQC (LB) | 13C2-4:2 Fluorotelomer sulfonate (13C2- 4:2 FTS) | 152 | 20 - 150 | 10 - 150 | percent | J/None | Y5 | No Qualifiers Applied |
| LABQC (LB) | 13C2-6:2 Fluorotelomer sulfonate (13C2- 6:2 FTS) | 156 | 20 - 150 | 10 - 150 | percent | J/None | Y5 | No Qualifiers Applied |

Where two qualifiers are listed, such as 'J/UJ', the first applies to positive results, and the second to non-detect results. Upper and Lower Warning and Control Limits are abbreviated UWL, LWL, UCL, and LCL in the Comment field.

No results associated with this QC element required qualification.

Data Validation Report for FC11160

Quality Control Outliers for test method E1633DR, LCS Recovery

The laboratory control sample/laboratory control sample duplicate (LCS/LCSD) serves as a monitor of the overall performance of each step during the analysis, including the sample preparation. Reported results were evaluated to determine compliance with the required acceptance criteria, and summary forms were evaluated and compared to electronic data deliverables. Findings of this review, and any associated qualified results, are listed below.

| Sample ID/ Lab Sample ID | Analyte | Result | Warning Limits | Control Limits | Units | Qualifier | Reason Code | Comment |
|-----------------------------|--|--------|-------------------|-------------------|---------|-----------|----------------|---------|
| OP58-BS (BS) | 3-Perfluoropropyl propanoic acid (3:3FTCA) | 180 | 40 - 150 | 40 - 150 | percent | J/None | C | |

Where two qualifiers are listed, such as 'J/UJ', the first applies to positive results, and the second to non-detect results. Upper and Lower Warning and Control Limits are abbreviated UWL, LWL, UCL, and LCL in the Comment field.

No results associated with this QC element required qualification.

Data Validation Report for FC11160

Qualified Results

No results associated with this sample delivery group required qualification.

Results with Modified Qualifiers

No qualifiers associated with this sample delivery group were modified manually.

Reason Code Definitions

| Code | Definition |
|------|-----------------------------|
| C | LCS Recovery |
| TR | Trace Level Detect |
| Y5 | Extracted Internal Standard |

Flag Code and Definitions

| Flag | Definition |
|------|---|
| J | Estimated Value |
| N | The analysis indicates the presence of an analyte for which there was presumptive evidence to make a tentative identification. |
| NJ | The analyte has been tentatively identified or presumptively as present and the associated numerical value was the estimated concentration in the sample. |
| R | The data are rejected due to deficiencies in meeting QC criteria and may not be used for decision making. |
| U | Undetected: The analyte was analyzed for, but not detected. |
| UJ | The analyte was not detected; however, the result is estimated due to discrepancies in meeting certain analyte-specific quality control criteria. |
| X | Result may require rejection; PDT attention required |

Bias

| | |
|---|-------------------------------|
| - | The result may be biased low |
| + | The result may be biased high |

Note - The bias field is a separate field; however, it is an integral part of the final flag (qualifier) on the sample result